

PRODUCT DATA SHEET

Standing: 2024-01-31


LAUDA PRO P 10 C

Heating bath thermostat 100-120 V; 50/60 Hz

Part Number: L000558

Features

- Thermostatic bath with latest microprocessor technology. Installed cooling coil
- Large 5,7" Multi-Touch-Color display
- Operation with multi touch display
- 360° turnable and detachable control panel for remote control
- User administration for up to 20 user
- Data in- and export via USB
- Timer und count-down-function
- EasyUse system for simple operation of the whole unit
- Low-level and adjustable over-temperature protection with acoustic alarm for use with flammable and non-flammable liquids
- SelfCheck Assistant for system diagnosis
- USB and Ethernet interface equipped as standard
- Retrofittable for upgrading with interface (RS 232/485, Profibus; EtherCAT; analogue or contact modules)
- Fully electronic continuous controller with PID action for internal & external control
- PowerAdapt system for the use of the maximum possible amount of heat permitted by the power supply system
- Programmer with up to 100 programs, up to 50 segments each
- Low height of appliance when control panel removed
- Significantly enlarged accessible volume
- Bath vessel made of stainless steel (insulated, with handles and drain tap)
- Internal LAUDA Varioflex pump with 8 selectable levels

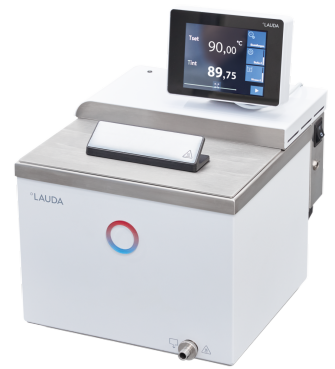

Working temperature min.
40 °C



Working temperature max.
250 °C

LAUDA DR. R. WOBSEY GMBH & CO. KG
Laudaplatz 1 • 97922 Lauda-Königshofen • DE

T + 49 (0) 9343 503-0
info@lauda.de • www.lauda.de
WEEE-Reg.-Nr.: DE 66 42 40 57

Kommanditgesellschaft: Sitz Lauda-Königshofen
Registergericht Mannheim • HRA 560069

Persönlich haftende Gesellschafterin:
LAUDA DR. R. WOBSEY Verwaltungs-GmbH
Sitz Lauda-Königshofen
Registergericht Mannheim • HRB 560226

Geschäftsführer:
Dr. Gunther Wobser (Vors.), Dr. Mario Englert,
Dr. Ralf Hermann, Dr. Marc Stricker
Beirat: Dr. Gerhard Wobser

PRODUCT DATA SHEET

Standing: 2024-01-31

LAUDA PRO P 10 C

Heating bath thermostat 100-120 V; 50/60 Hz

Part Number: L000558

Technical Features (according to DIN 12876)

Working temperature range	40 ... 250 °C
Working temperature range with external cooling	20 ... 250 °C
Operating temperature range	-30 ... 250 °C
Ambient temperature range	5 ... 40 °C
Temperature stability	0.01 ± K
Heater power max.	1.85 kW
Power consumption max.	1.9 kW
Current max.	16 A
Size of bath (W x D x H)	240 x 150 x 200 mm
Bath volume min. / max.	5.0 / 10.0 L
Overall dimensions (WxDxH)	310 x 335 x 415 mm
Power supply	100-120 V; 50/60 Hz
Power plug	Power cord without plug (AWG)

Standard accessories

- 1 Bath cover
- 2 nipples 10 mm with screw cap G3/8 for cooling water

LAUDA DR. R. WOBSEY GMBH & CO. KG
Laudaplatz 1 • 97922 Lauda-Königshofen • DE

T + 49 (0) 9343 503-0
info@lauda.de • www.lauda.de
WEEE-Reg.-Nr.: DE 66 42 40 57

Kommanditgesellschaft: Sitz Lauda-Königshofen
Registergericht Mannheim • HRA 560069

Persönlich haftende Gesellschafterin:
LAUDA DR. R. WOBSEY Verwaltungs-GmbH
Sitz Lauda-Königshofen
Registergericht Mannheim • HRB 560226

Geschäftsführer:
Dr. Gunther Wobser (Vors.), Dr. Mario Englert,
Dr. Ralf Hermann, Dr. Marc Stricker
Beirat: Dr. Gerhard Wobser


PRODUCT DATA SHEET

Standing: 2024-01-31

LAUDA PRO P 10 C

Heating bath thermostat 100-120 V; 50/60 Hz

Part Number: L000558

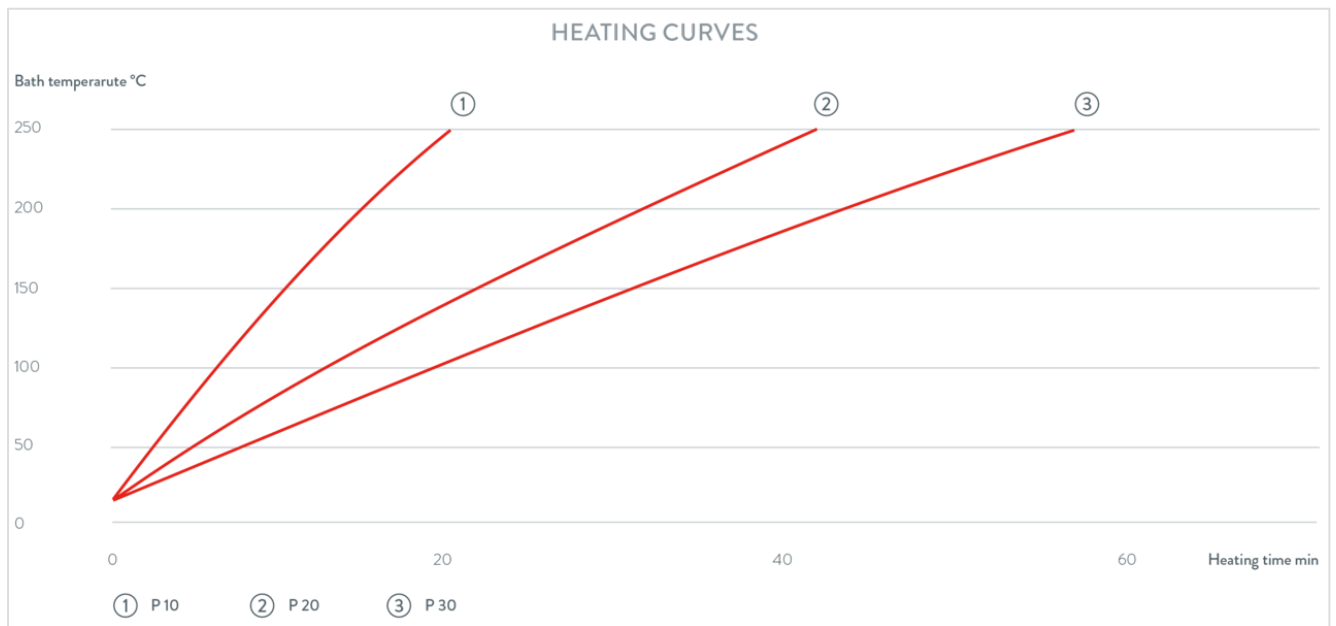

Reserve technical changes

LAUDA DR. R. WOBSEY GMBH & CO. KG
Laudaplatz 1 • 97922 Lauda-Königshofen • DE

T + 49 (0) 9343 503-0
info@lauda.de • www.lauda.de
WEEE-Reg-Nr.: DE 66 42 40 57

Kommanditgesellschaft: Sitz Lauda-Königshofen
Registergericht Mannheim • HRA 560069

Persönlich haftende Gesellschafterin:
LAUDA DR. R. WOBSEY Verwaltungs-GmbH
Sitz Lauda-Königshofen
Registergericht Mannheim • HRB 560226

Geschäftsführer:
Dr. Gunther Wobser (Vors.), Dr. Mario Englert,
Dr. Ralf Hermann, Dr. Marc Stricker
Beirat: Dr. Gerhard Wobser